## TEMPERATURE AND SPEED CONTROL LAB (TSC-LAB)

Práctica 5: Activation of Transistor 1 and 2, also Reading of temperature sensor 1 and 2

https://tsc-lab.blogspot.com/2021/05/practica-5-activation-of-transistor-1.html

```
    Blog: https://tsc-lab.blogspot.com/
```

- GitHub: https://github.com/vasanza/TSC-Lab
- Matlab functions: https://github.com/vasanza/Matlab Code
- IEEEDataPort: http://ieee-dataport.org/4138
- TSC-LAB configurations
- int period = 15; //medium period in minutes
- int freq\_sampling = 100; // sampling time
- int ciclos = 20; // sampling time
- 1 sample per second

## Raw dataset preparation


```
clear;clc;%clear all
addpath(genpath('./src'))%functions folders
datapath = fullfile('./data/');%data folder
```

### Raw dataset preprocessing

```
filenames = FindCSV(datapath);%List All CSV files
data=readtable(fullfile(datapath,filenames(1).name));%Select i CSV file
data=table2array(data);
DataNorm = fNormalization(data(:,1:2));%Normalization
DataFeatures = [max(DataNorm) min(DataNorm) mean(DataNorm)...
    median(DataNorm) rms(DataNorm) std(DataNorm)];%Feature extraction
```


#### Plot Raw TSC-LAB dataset

```
figure
plot(data);xlabel('Samples');ylabel('°C');
title('Temperature');
legend('Temp1','Temp2','PWM1','PWM2');
```

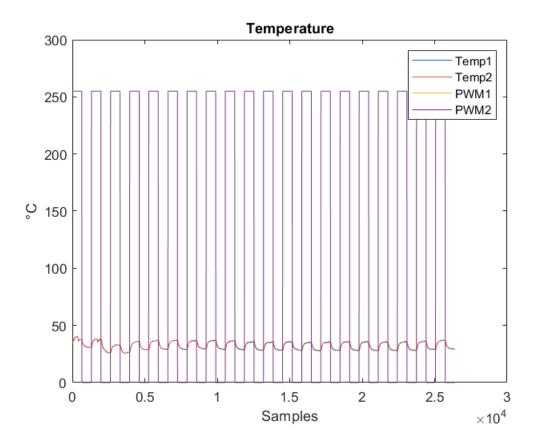


# Plot Normalization dataset

```
figure
plot(DataNorm);xlabel('Samples');ylabel('Normalized Values');
title('Normalized data graph');
legend('Temp1','Temp2');
```

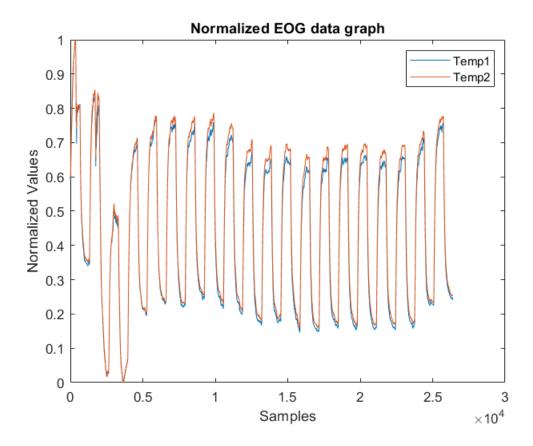

# Select a case

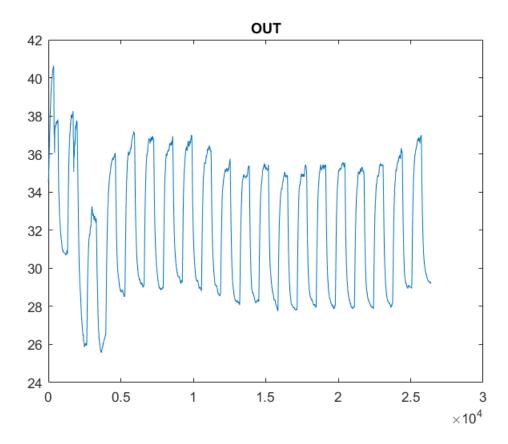


# **System Identification**

#### ident

Warning: The "ident" command is obsolete and may be removed in a future release of MATLAB. Use the "systemIdentification" command instead.

Created preference file C:\Users\vasan\Documents\MATLAB\idprefs.mat. Type HELP MIDPREFS if you want to move this file.

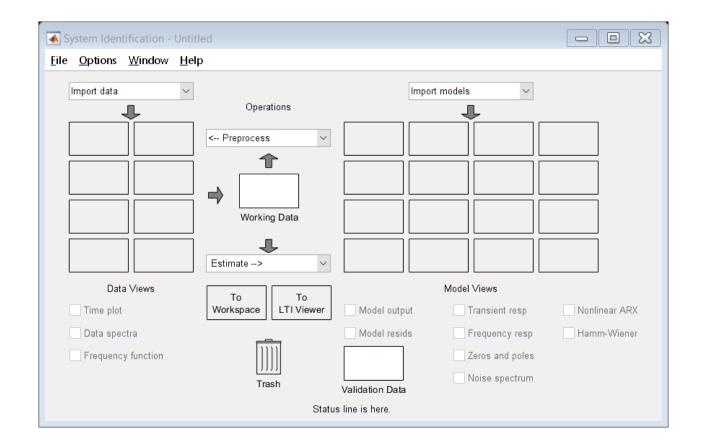

# Open the Classification Learner

%regressionLearner
%classificationLearner